Orgo-Life® Rife Frequency Generator UL-9051



Brand: Orgo-Life® Technology
Product Code: UL-9051
Availability: 5
Weight: 0.00kg
Dimensions: 0.00cm x 0.00cm x 0.00cm

Price: $189.95

Ex Tax: $189.95

Short Description

Orgo-Life® Rife Frequency Generator UL-9051—Simple , fast & Ideal for beginners doing Bio-Resonance

The **Orgo-Life® UL-9051 Frequency Generator** is an enhanced, stabilized version designed specifically for the Orgo-Life® ecosystem. It works in full harmony with **all vibrational units, quantum modules, coils, antennas, energetic plates, and wellness test systems**.

Description

UL-9051 Rife System – 1 Hz to 500 kHz Frequency Generator + UL7013 Interface + Full Orgo-Life® Accessories

quantum modules, coils, antennas, energetic plates, and wellness test systems.

Everything is **100% non-medical**, intended **exclusively for wellness purposes**, fully compliant with **Law 31**.

This system allows you to perform:

- personal tests
  - harmonization frequency sessions
  - wellness-exploration tests

• vibrational experiences to restore inner balance
• reduction of “bad vibes”
• calming internal tensions
• regaining energetic stability in the body

The generator is officially delivered with **all Orgo-Life® cables and interfaces**, ready to connect to every UL device, spiral, module, plate, coil and Orgo-Life® accessory.

1 Official Orgo-Life® Kit Content

Each system includes:

— Orgo-Life® UL-7013 Interface

Allows connection to:
• UL modules
• spirals
• headsets
• plates
• discs
• antennas
• copper tubes
• Orgo-Life® quantum accessories

— Mini 1/8 Interface (included)

A compact universal adapter allowing connection to **all Orgo-Life® units**, including older generations.

— Orgo-Life® Tesla Coil Antenna (your choice)

- 50 mm
  • 100 mm
  • 150 mm

Each coil includes:
? dedicated cable
? Orgo-Life® adapter
? UL-7013 interface connection support

— Orgo-Life® Copper Tubes (2 units)

- Length: 6 inches
  • Diameter: ¾ inch
  • Material: 100% pure copper
  Used for hand-held frequency tests or harmonization sessions.

— Orgo-Life® BNC Alligator Clip Cables

Used to connect:
• UL modules
• plates
• copper tubes
• antennas
• various vibrational test setups

— Additional Coil Cable (50 / 100 / 150 mm depending on your choice)

Designed to transmit frequencies with stability.

— High-stability USB power cable (3C certified)

2 General Description of the UL-9051 Orgo-Life® Generator

This DDS (Direct Digital Synthesis) generator can produce **over 15 vibrational waveform types**, all compatible with Orgo-Life® systems.

Available Waveforms:

- Sine
  • Square
  • Triangle
  • Sawtooth + Reverse Sawtooth
  • ECG
  • Schumann (7.83 Hz default)
  • Gaussian
  • Lorentz
  • Sinc
  • Exponential rise
  • Exponential decay
  • Half-wave
  • Full-wave
  • Noise
  …and all standard international waveform variations used in vibrational

testing.

Technical Specifications:

- **Frequency range:**
  – Sine: 0.01 Hz ? 500 kHz
  – Other waveforms: 0.01 Hz ? 100 kHz
  • **Resolution:** 0.01 Hz (ultra-precise)
  • **Max amplitude:** ±10V (Peak-to-Peak)
  • **Impedance:** 50 ohms
  • **DC Offset:** ±10V (switchable)
  • **Power supply:** USB DC 5V

Integrated Filters

Designed to produce cleaner signals, especially for:
• high-frequency sine waves
• pulsed signals
• Rodin coils / Orgo-Life® antennas

Typical Orgo-Life® Applications:

- Harmonization
  • Vibrational cleansing
  • Personal quantum exploration
  • Emotional stabilization (wellness only, non-medical)
  • Subtle meridian realignment
  • Mental calming
  • Exploratory frequency testing

3 Controls and Functions

MODE — Changes the waveform
CURSOR — Selects the digit to edit
+ / ? — Adjusts the frequency
RUN/STOP — Starts or stops output
DCOFFSET — Activates/deactivates DC component
FILTER — Cleans the output signal
AMP — Adjusts signal amplitude

4 Steps for Using With Orgo-Life®

1. Plug the USB power cable into the generator.
2. Connect the BNC cable to the UL-7013 interface.
3. Connect your modules:
   coil antenna
   • pure copper tubes (¾ inch × 6 inches)
   • Tesla spiral 50 mm, 100 mm or 150 mm
   • quantum plate
4. Select your waveform.
5. Enter a frequency (e.g., 7.83 Hz, 432 Hz, 963 Hz).
6. Press **RUN**.

5 Wellness Application (Law 31 Compliant)

This generator is used for:
• centering oneself
• reducing negative energy or internal discomfort
• regaining a stable state
• supporting deep relaxation
• improving vibrational sensitivity
• exploring frequencies that help “clear the mind”
• harmonizing intentions and emotions

Always for **wellness purposes only**, never medical, in accordance with **Law 31**.

6 Mandatory Legal Notice (Law 31)

Orgo-Life® devices are intended exclusively for wellness, harmonization and energetic relaxation.
They **do not diagnose, treat, or prevent any disease**.
They do not replace any healthcare professional.
For any physical or psychological condition, consult a licensed medical specialist.

TECHNICAL SPECIFICATIONS

MODEL: UL-9051

DESCRIPTION:
• 3 × downloadable PDF booklets (on Orgo-Life.com) with product description (FR/EN)
• 1 × USB power block
• 1 × USB charging cable
• 1 × alligator cable (optional) with 2 × ¾ inch pure copper pipes

• 1 × UL-7013 interface cable to plug Orgo-Life® units
• 1 × Tesla antenna 50 mm, 100 mm or 150 mm with cable and adapter plug
• 1 × BNC to mini 1/8 cable for FI headphone (very comfortable)

VOLTAGE:
110 volts or 220 volts – all devices are universal.
All devices come with 220 V adapters (for countries outside North America), compatible with French/European current, and work on 240 V as well as 110 V.

PACKAGING:
35 cm × 25 cm × 15 cm – approx. 3 kg

BARCODE:
787416647851

INSTRUCTIONS:
Comes with French & English instruction manuals in PDF (downloadable on Orgo-Life.com)

DISCLAIMER:
Our devices do not claim to cure any disease or illness.
In the event of a serious medical condition, please seek the help and advice of a qualified medical professional.
The potential health benefits have not been evaluated by national medical licensing bodies.

MANUFACTURER:
Original Health Distribution Ltd., Canada
Email: original.health.net@gmail.com

PRODUCT SAFETY:
Instructions are downloadable on the Orgo-Life.com website (PDF document).

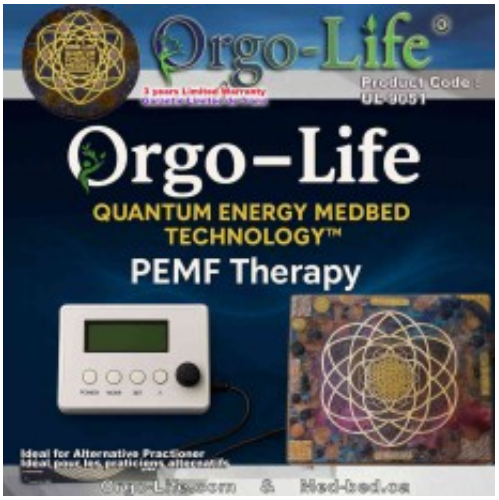
PRODUCTION PARTNERS:
OrgoLifeTechnology, with support from Renewable Free Energy Foundation, Magog, QC, Canada.

Product Gallery


Orgo-Life
Product Code :
UL-9051
1 USB Transformer
1-DC connector cable
1 Tesla Antenna 50 mm + Adaptor mini 1/8 plug
1 X UL-7013 cable for Orgo-Life Unit
-1 cable for Antenna mini 1/8 , 1 BNC to Female 1/8
-1 Cable to bare alligator plugs for Copper pipes
2 X copper pipe 3/4 inch diameter X 6 inches approx.
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-9051
Ideal for Alternative Practioner
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-9051
Ideal for Alternative Practioner


Orgo-Life
Product Code :
UL-9051
Orgo-Life
QUANTUM ENERGY MEDBED TECHNOLOGY™
PEMF Therapy
Ideal for Alternative Practioner
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-9051
Orgo-Life
QUANTUM ENERGY
ENERGY TECHNOLOGY®
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-9051
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-9051
QUANTUM
ENERGY
MEDBED
TECHNOLOGY®
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.us


Orgo-Life
Product Code :
UL-9051
Orgo-Life.com & Med-bed.us


Orgo-Life
Product Code :
UL-9051
Use the copper pipes options to do special
Rife Frequency work on your body
Orgo-Life.com & Med-bed.us


Orgo-Life
Product Code :
UL-9051
Compatible with all Orgo-life® Units
Also compatible with Quantum Boosters
UL-8055 (1000%) & UL-8050 (300%)
Orgo-Life.com & Med-bed.us


Orgo-Life
Product Code :
UL-9051
Option without copper pipes & its cable
With Tesla Antenna 50 mm + cables
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs


Orgo-Life
Product Code :
UL-9051
Option without copper pipes & its cable
With Tesla Antenna 100 mm + cables
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs


Orgo-Life®
Product Code :
UL-9051
Option without copper pipes & its cable
With Tesla Antenna 150 mm + cables
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs


Orgo-Life®
Product Code :
UL-9051
1 USB Transformer
1-DC connector cable
1 Tesla Antenna 100 mm + Adaptor mini 1/8 plug
1 X UL-7013 cable for Orgo-Life Unit
-1 cable for Antenna mini 1/8 , 1 BNC to Female 1/8
-1 Cable to bare alligator plugs for Copper pipes
-2 X copper pipe 3/4 inch diameter X 6 inches approx.
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life®
Product Code :
UL-9051
1 USB Transformer
1-DC connector cable
1 Tesla Antenna 150 mm + Adaptor mini 1/8 plug
1 X UL-7013 cable for Orgo-Life Unit
-1 cable for Antenna mini 1/8 , 1 BNC to Female 1/8
-1 Cable to bare alligator plugs for Copper pipes
-2 X copper pipe 3/4 inch diameter X 6 inches approx.
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life®
Product Code :
UL-9051
QUANTUM
ENERGY
MEDBED
TECHNOLOGY®
Ideal for Alternative Practioner
Ideal pour les praticiens alternatifs
Orgo-Life.com & Med-bed.ca


Orgo-Life®
Product Code :
UL-9000
UL-7013 interface cable to
plug any compatible
Orgo-Life Unit
Size: 150MM
Available in 150mm 100mm
or 50mm Tesla Antenna
Orgo-Life.com & Med-bed.ca


Orgo-Life®
Product Code :
UL-9051
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-8051
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-8051
Orgo-Life
QUANTUM ENERGY
ENERGY TECHNOLOGY®
Orgo-Life.com & Med-bed.ca


Orgo-Life
Product Code :
UL-8051
Square
Triangle
100mm/3.94"
150mm/5.12"


Orgo-Life
Product Code :
UL-8051
Square
Triangle
100mm/3.94"
150mm/5.12"


Orgo-Life
Product Code :
UL-8051
Compatible with all Orgo-life® Units
Also compatible with Quantum Boosters
UL-8055 (1000%) & UL-8050 (300%)
UL-8030
Orgo-Life.com & Med-bed.ca